



<http://www.unisis.co.jp>

UNIEVER[®]

Uniever is a brand name of Unisis Corp.


DISPOSABLE NERVE BLOCKADE NEEDLE



25G×25mm
(for SGB)



22G×60mm

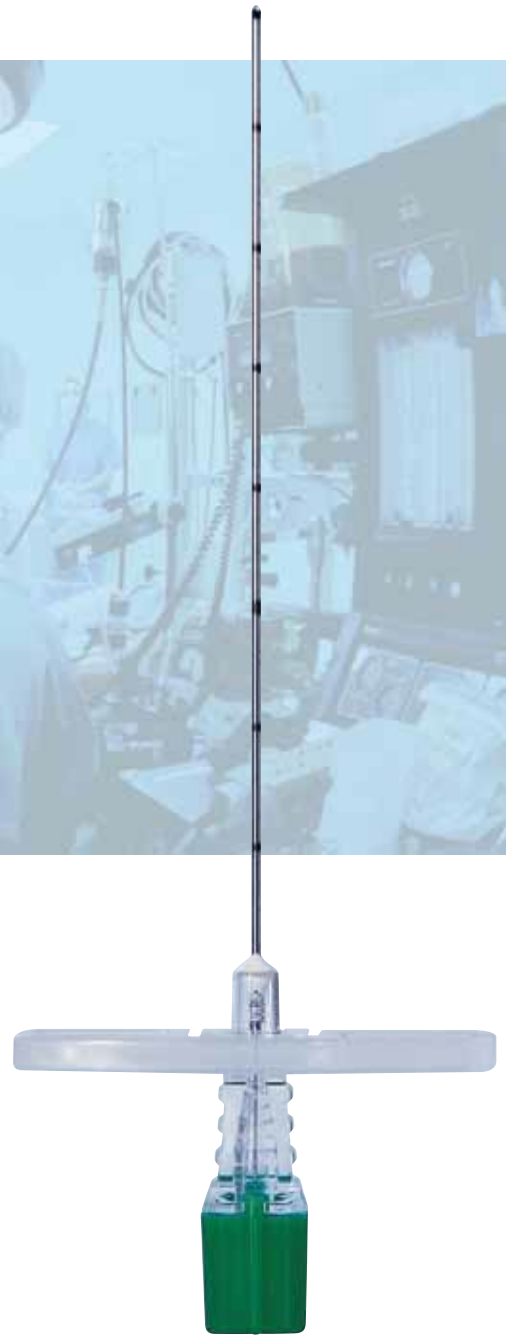

21G×80mm

DISPOSABLE NERVE BLOCKADE NEEDLE

Indispensable needles for pain Clinics, based on the real needs of physicians.

Our nerve blockade needles were developed to provide the cutting performance that all physicians need. We have two different point types that focus on puncture site and detection. A wide variety of size options are available for blocking different nerve ganglions and plexuses.

Point Symbol Marks



Back-Cut point - short bevel



K-3 Lancet point - short bevel

Features

● Ideal point shapes for the purpose

Nerve blockade needles should provide two adverse features, smooth puncture and appropriate resistance. Unisis provides two different types of point shape to attain these goals, namely the K-3 Lancet point and the Back-Cut point, so that the surgeon can choose based on the requirements of the puncture site.



● Precise fitting of stylet and needle

The clearance between stylet and needle is consistently minimized by our strict quality control. The precise bevel facet formation reduces damage to the tissue and its cutting.

● Crystal-clear hub

The crystal-clear polycarbonate hub allows quick recognition of an unexpected hemorrhage. The wing is detachable to support various uses.

22G needles w/different points puncture patterns with nylon film.

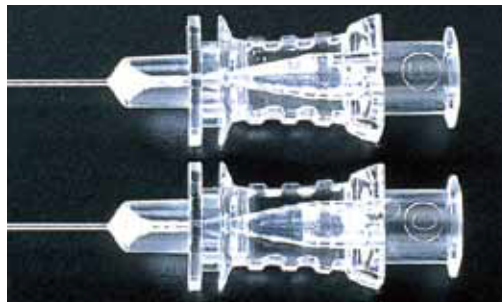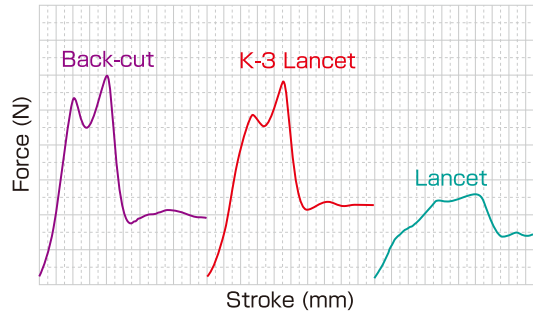

Needle points



Back-Cut point (short bevel)


Standard used point in nerve block needles, which is characterized by an excellent straight puncture and gives an adequate feel of the nerve plexus. A short-beveled secondary blade (bevel angle: 30 degrees) is cut on each side to give an adequate feel for recognizing insertion into the nerve plexus.



K-3 Lancet point (short bevel)

Our short-beveled point (bevel angle: 30 degrees) supports safe puncture. With cutting performance superior to that of the Back-Cut point, this short-beveled needle can be used in the same way as ordinary Lancet point needles.

Names of parts with sketches

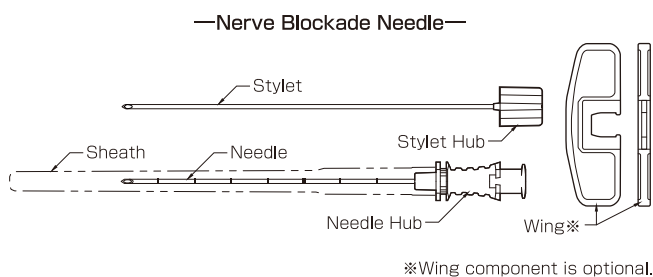

Standard sizes

Product Name : UNIEVER DISPOSABLE NERVE BLOCKADE NEEDLE

Needle Point Type : Back-Cut Point

G: Gauge/ Color Code	Point type	Length	Graduation	Attachment	REF (Code No.)
22G	Back-Cut	60mm	Line Type	Wing	4370771
		80mm	Line Type	Wing	4371471
		100mm	Line Type	Wing	4372071

Individually packed and 25 pcs/box
E.O.G. sterilized

Manufacturer



7th Fl. Mitsui-Sumitomo Bank Okachimachi-Bldg. 4-11-4 Taito, Taito-Ku,
Tokyo, Japan 110-0016
(Direct call) International Sales Dept. Tel: +81-3-5812-8832 / Fax: +81-3-5812-8831

URL <http://www.unisis.co.jp/english/>

Made in Japan

Representative

Note : Product specifications are subject to change without prior notice due to product improvements, etc. NB090104E