



<http://www.unisis.co.jp>

UNIEVER[®]

Uniever is a brand name of Unisis Corp.

DISPOSABLE EPIDURAL ANESTHESIA MINITRAY (Pediatric)



DISPOSABLE EPIDURAL ANESTHESIA MINITRAY FOR PEDIATRIC USE



Huber point

The excellent entry kit for epidural anesthesia designed to fit pediatric body.

Devices required for pain management during and after operation are assembled in a compact kit. Combination of flexible easy-to-handle catheter and epidural needle fitting to pediatric patient forms can support your dedicated procedure during pediatric epidural anesthesia.

Characteristics

●Catheter for Epidural Anesthesia


〈Biocompatibility〉

22G Polyurethane catheter that fits to pediatric body has an excellent biocompatibility to have temperature sensitive hardness, being hard enough at insertion into body; it will be softened by body temperature and minimize risk of breakage.



〈Easy detection (Radiopaque line)〉

Radiopaque agent is linearly embedded for easy detection of catheter piece left, in case of breakage inside body.

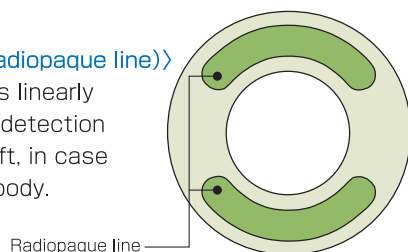


Radiopaque line

●Epidural Anesthesia Needle (Huber)

〈Anti-coring process at bevel-heel〉

Heel part of the bevel is specially treated for preventing catheter damage such as breakage.

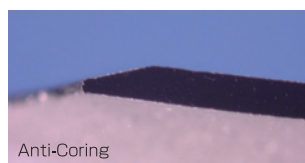

Anti-Coring



Without Anti-Coring

〈Epidural anesthesia needle/Tip markings〉


The standard sizes of epidural anesthesia needle include 19G/35mm and 19G/50mm, which fit to pediatric body. Upon requests of physicians giving anesthetic into pediatric patients, the needle has depth markings. This allows you to carefully manipulate the epidural needle.




Kit Components

●Catheter Connector

Ergonomically designed connector offers you firm holding and easy handling of catheter provided by screwing mechanism. Our unique technology delivered the inner structure of catheter, much less likely to cause obturation (patent issued).

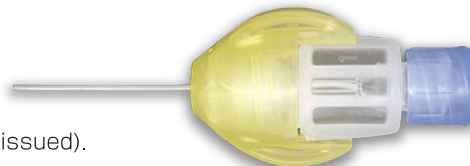

Kit components (basic)

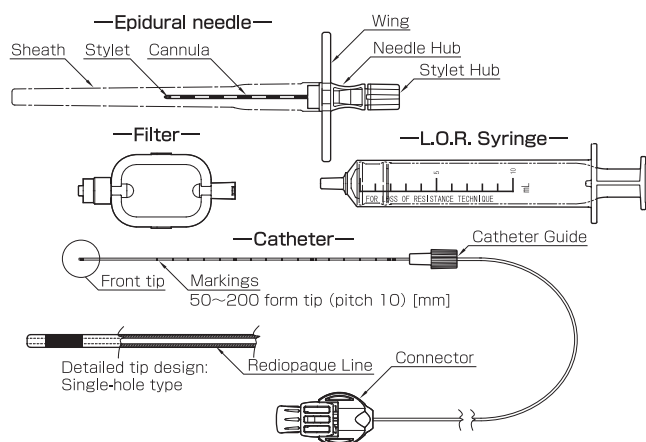

Standard kit

Product Name : DISPOSABLE EPIDURAL ANESTHESIA MINITRAY (Pediatric)

Epidural Anesthesia Needle Huber Point	Epidural Catheter	L.O.R. Syringe	Accessory	Product code	
G:Gauge/Color Code	Length				
19G	35mm	Single-hole type	10mL	·Catheter Guide ·Catheter Connector ·Filter	5143261
	50mm			·Epidural label One per piece	5143361

Individually packed and 10 pcs/box
E.O.G. sterilized

*For more information on products other than the standard products listed above, please contact the sales department in charge.

Manufacturer

UNISIS CORP. JAPAN

7th Fl. Mitsui-Sumitomo Bank Okachimachi-Bldg. 4-11-4 Taito, Taito-Ku,
Tokyo, Japan 110-0016
(Direct call) international sales dept Tel:03-5812-8832 / Fax:03-5812-8831

URL <http://www.unisis.co.jp/english/>

Made in Japan

Representative