



<http://www.unisis.co.jp>

UNIEVER[®]

Uniever is a brand name of Unisis Corp.

DISPOSABLE SPINAL ANESTHESIA NEEDLE (K-3 LANCET POINT)



29G×75mm

26G×90mm

24G×120mm

DISPOSABLE SPINAL ANESTHESIA NEEDLE (K-3 LANCET POINT)

Point Symbol Marks



K-3 Lancet point



K-3 Lancet point
(with Introducer Needle)

Spinal needle continues to set the global standard in a rapidly evolving category.

This is standard spinal needle with conventional Quincke point, featuring excellent cutting performance and easy manipulation, ensuring successful procedures. It is available in a wide variety of sizes at reasonable prices to support different purposes and procedures.


Features

● Crystal-clear hub and rapid flush back of CSF

Our crystal-clear polycarbonate hub allows easy visualization of CSF flush back. Scientifically designed hub provides a secure grip during procedure.

A carefully designed and molded hub provides world class performance for the rapid flush back of CSF.


Our strict QC procedures guarantee perfect fitting of stylet and needle hubs, which prevents needle wobbling.

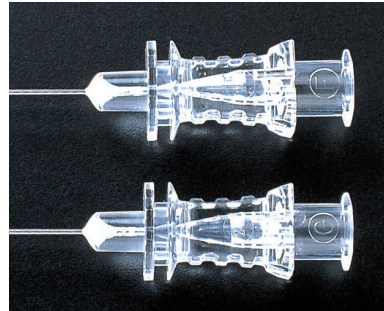

● Bevel coordination


Slight deviation of stylet/needle bevel angle is intended to minimize dural damages.

We pay careful attention to the assembly of the needle and stylet in order to avoid problems like tissues jammed between the bevels.

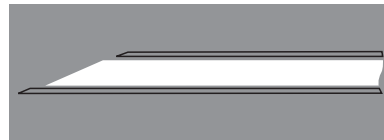


● Thin wall

Our recent application of unique thin-walled materials has resulted in an inner diameter 30% larger compared to our previous model. You'll see consistently fast CSF flush back thanks to this maximum lumen and our innovative hub design.

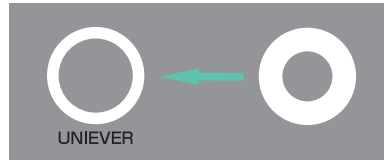

Needle points

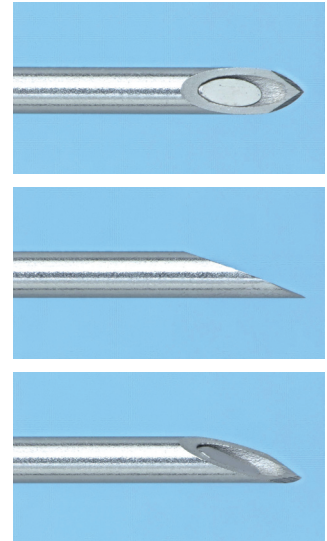


K-3 Lancet point


Our K-3 Lancet (modified Quincke) point with secondary side bevels ground at appropriate angle provides moderate puncture resistance, while minimizing dural damages.

Names of parts with sketches

—Spinal Anesthesia Needle (K-3 Lancet Point)—

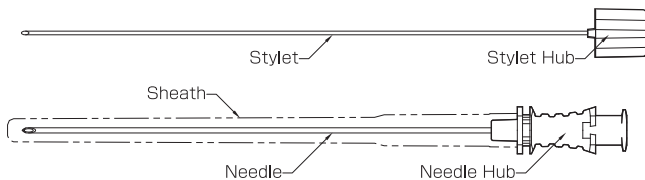

—Introducer—

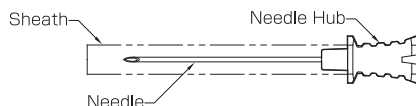

Standard sizes

Product Name : UNIEVER DISPOSABLE SPINAL ANESTHESIA NEEDLE
Needle Point Type : K-3 Lancet Point

G:Gauge/ Color Code	Length	Introducer	REF(Code No.)
20G	90mm	—	1131161
22G	90mm	—	1151261
23G	90mm	—	1160961
25G	90mm	—	1182961
26G	90mm	20G38mm	5240561
27G	90mm	22G38mm	5250561
29G	90mm	22G38mm	5260461

Individually packed and 25 pcs/box
E.O.G. sterilized

Manufacturer



7th Fl. Mitsui-Sumitomo Bank Okachimachi-Bldg. 4-11-4 Taito, Taito-Ku,
Tokyo, Japan 110-0016
(Direct call) International Sales Dept. Tel: +81-3-5812-8832 / Fax: +81-3-5812-8831

URL <http://www.unisis.co.jp/english/>

Made in Japan

Representative

Note : Product specifications are subject to change without prior notice due to product improvements, etc. **SP090204E**